

Anti - Ki-67

Rabbit clonal antibody

CAT#

CONCENTRATED

DB 070-0.1	(100 µl)
DB 070-0.2	(200 µl)
DB 070-0.5	(500 µl)
DB 070-1	(1 ml)

READY TO USE (RTU)

DB 070-RTU-7	(7 ml)
DB 070-RTU-15	(15 ml)

STORAGE AND APPLICATION

CONCENTRATED

Storage: +4°C
Application: IHC-P,
dilution 1:100 - 1:200

READY TO USE (RTU)

Storage: +4°C, Do not freeze!
Application: IHC-P,
ready to use

PRODUCT INFORMATION

Clone: C16-I
Buffer: 20 mM Tris-HCl, pH 8.0
Stabilizer: 20 mg/ml BSA
Preservative: 0.05% NaN₃
Specificity: Human
Expiration: 24 months from the shipping date
Immunogen: Peptide derived from C-terminal region of human Ki-67. Antibody recognizes the epitope between Lys3243 - Asp3255.

Cellular localization: nucleus
Positive control: tonsil
Protein accession number: P46013


IHC-P PROTOCOL – INSTRUCTION MANUAL

1. Deparaffinize the section in 3 changes of xylene, 5 minutes each.
2. Wash the section in 96%, 80% and 70% ethyl alcohol for 10 minutes each.
3. Rinse in distilled water.
4. Block the endogenous peroxidase by incubating the tissue in 3% hydrogen peroxide (H₂O₂) for 10 minutes.
5. Wash in distilled water.
6. For antigen retrieval: immerse the slide in Tris-EDTA buffer, pH 9.0, 0.05% Tween-20*, and incubate in water bath at 96-98°C for 20-25 minutes. (Alternatively adjust to your own protocol, keeping the required pH)
7. Transfer the slide to room temperature and let it cool down (in Tris-EDTA buffer, pH 9.0) for 15 minutes.
8. Rinse in distilled water.
9. Wash in 0.05 M Tris-HCl, pH 7.6 buffer supplemented with 0.2% of Tween-20 (**buffer A**) for 5 minutes.
10. **CONCENTRATED:**
Incubate the section with primary antibody at the **dilution 1:100 - 1:200** for 1 hour in the closed wet chamber.
READY TO USE (RTU):
Incubate the section with primary antibody (**ready to use**) for 1 hour in a closed wet chamber.
11. Wash twice 5 minutes with buffer A.
12. Apply the secondary antibody (the protocol depends on the supplier), and proceed to standard immunohistochemistry protocol (HRP - Peroxide - DAB).
Micropolymer-HRP detection kit rabbit/mouse dual of DB Biotech is suggested (<http://www.dbbiotech.com/products/detection-system.html>).
13. Wash twice 5 minutes with buffer A.
14. Apply the chromogen (DAB), 1 - 3 minutes.
15. Wash in water – 10 minutes.
16. Stain in hematoxylin for 5 minutes.
17. Wash in water – 10 minutes.
18. Dehydrate the section in 2 changes of 96% ethyl alcohol for 5 minutes each.
19. Wash the section in 2 changes of xylene for 2 minutes each.
20. Mount the slide for observation.

*Tris-EDTA Buffer (10mM Tris Base, 1mM EDTA solution, 0.05% Tween-20, pH 9.0):
Tris ----- 1.21 g; EDTA ----- 0.37 g; Distilled water ----- 1000 ml
Mix to dissolve in 700 ml of distilled water. Adjust pH to 9.0 with 1M HCl and then add 0.5 ml of Tween-20 and mix well. Adjust the final volume to 1 liter with distilled water.

PRECAUTIONS

1. **We strongly recommend to use DB Primary Antibody Diluent (catalog number DB D-125, or DB D-250), eventually alternative diluent (containing protease free BSA at the concentrations ≥ 1mg/ml) for dilution of concentrated antibodies, otherwise the warranty might be voided.**
2. **Centrifuge the vial before use.**
3. Intended for professional In Vitro Diagnostic use in laboratories.
4. Do not use after expiration date stamped on vial label.
5. Avoid contamination of the reagent.
6. Any discrepancies in the recommended procedures stated in the working protocol may affect the final results.
7. The reagent contains sodium azide (NaN₃) which is highly toxic in higher concentrations. The concentration in the reagent (0.05%) is not considered as hazardous.
8. Disposal of waste material must be conducted in accordance with local regulations.
9. Wear appropriate Personal Protective Equipment to avoid contact with eyes and skin.

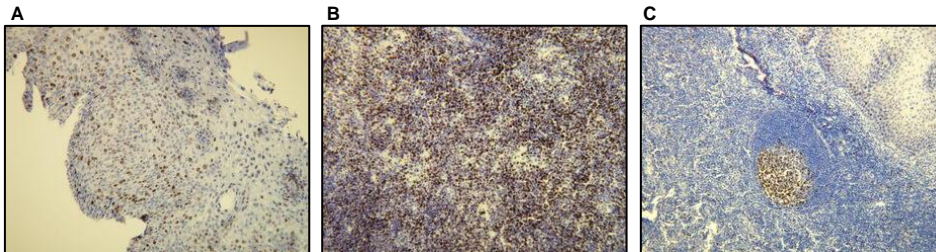

High Grade SIL (squamous intraepithelial lesion) of the uterine cervix showing full-thickness Ki-67 immunopositivity in the dysplastic epithelium (A), B-cell lymphoma showing high proliferating index visualised by the Ki-67 immunohistochemistry (B) and proliferating germinal centre cells of the lymphoid follicle visualised by the Ki-67 immunohistochemistry (C). All, formalin fixed, paraffin embedded human tissues (4 µm sections) stained with anti - Ki-67 (DB 070) monospecific clonal antibody according to related DB Biotech datasheet.