

Anti - Alpha Smooth Muscle Actin

Rabbit clonal antibody

CAT#

CONCENTRATED

DB 147-0.1	(100 µl)
DB 147-0.2	(200 µl)
DB 147-0.5	(500 µl)
DB 147-1	(1 ml)

READY TO USE (RTU)

DB 147-RTU-7	(7 ml)
DB 147-RTU-15	(15 ml)

STORAGE AND APPLICATION

CONCENTRATED

Storage: +4°C
Application: IHC-P,
 dilution 1:100 - 1:200

READY TO USE (RTU)

Storage: +4°C, Do not freeze!
Application: IHC-P,
 ready to use

PRODUCT INFORMATION

Clone: M16-L
Buffer: 20 mM Tris-HCl, pH 8.0
Stabilizer: 20 mg/ml BSA
Preservative: 0.05% NaN₃


Specificity: Human, mouse
Expiration: 24 months from the shipping date
Immunogen: Peptide derived from the N-terminal sequence of human alpha smooth muscle actin. Antibody recognizes the epitope between Glu3 - Gly16.

Cellular localization: cytoplasm, cytoskeleton
Positive control: tumors arising from smooth muscles and myoepithelial cells
Protein accession number: P62736

IHC-P PROTOCOL – INSTRUCTION MANUAL

1. Deparaffinize the section in 3 changes of xylene, 10 minutes each.
2. Wash the section in 96%, 80% and 70% ethyl alcohol for 10 minutes each.
3. Rinse in distilled water, 2 x 5 minutes.
4. Block the endogenous peroxidase by incubating the tissue in 3% hydrogen peroxide (H₂O₂) for 10 minutes.
5. Wash in distilled water, 2 x 5 minutes.
6. For antigen retrieval: Immerse the slide in Tris-EDTA buffer, pH 9.0* and incubate at 95-97°C in water bath for 20 minutes.
7. Remove the staining to room temperature and let the slide to cool (in Tris-EDTA buffer, pH 9.0) for 15 minutes.
8. Rinse in distilled water, 2 x 5 minutes.
9. Wash in PBS (phosphate buffer saline, pH 7.0-7.5) supplemented with 0.05% of Tween-20 (**buffer A**), 2 x 5 minutes
10. **CONCENTRATED:**
 Incubate the section with primary antibody at the **dilution 1:100 – 1:200** for 1 hour in the closed wet chamber.
READY TO USE (RTU):
 Incubate the section with primary antibody (**ready to use**) for 1 hour in a closed wet chamber.
11. Wash 3 x 5 minutes with buffer A.
12. Apply the secondary antibody (the protocol depends on the supplier), and proceed to standard immunohistochemistry protocol (HRP - Peroxide - DAB). Micropolymer-HRP detection kit rabbit/mouse dual of DB Biotech is suggested (<http://www.dbbiotech.com/products/detection-system.html>).
13. Wash 3 x 5 minutes with buffer A.
14. Apply the chromogen (DAB), 1 - 3 minutes.
15. Wash in water, 2 x 5 minutes.
16. Stain in hematoxylin for 5 minutes.
17. Wash in distilled water, 3 x 2 minutes.
18. Mount the slide for observation.

* **Tris-EDTA Buffer (10mM Tris Base, 1mM EDTA solution, 0.05% Tween-20, pH 9.0):**
 Tris ----- 1.21 g; EDTA ----- 0.37 g; Distilled water ----- 1000 ml
 Mix to dissolve in 700 ml of distilled water. Adjust pH to 9.0 with 1M HCl and then add 0.5 ml of Tween-20 and mix well. Adjust the final volume to 1 liter with distilled water.
 Store this solution at room temperature for 3 months or at +4°C for longer storage.

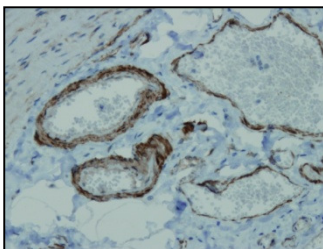

Alpha smooth muscle actin positivity in the vessel walls of submucosa of small intestine tissue. Formalin fixed, paraffin embedded human tissue (4 µm section) stained with DB 147 monospecific antibody according to related DB Biotech datasheet.

PRECAUTIONS

1. **We strongly recommend to use DB Primary Antibody Diluent (catalog number DB D-125, or DB D-250), eventually alternative diluent (containing protease free BSA at the concentrations ≥ 1mg/ml) for dilution of concentrated antibodies, otherwise the warranty might be voided.**
2. **Centrifuge the vial before use.**
3. Intended for professional In Vitro Diagnostic use in laboratories.
4. Do not use after expiration date stamped on vial label.
5. Avoid contamination of the reagent.
6. Any discrepancies in the recommended procedures stated in the working protocol may affect the final results.
7. The reagent contains sodium azide (NaN₃) which is highly toxic in higher concentrations. The concentration in the reagent (0.05%) is not considered as hazardous.
8. Disposal of waste material must be conducted in accordance with local regulations.
9. Wear appropriate Personal Protective Equipment to avoid contact with eyes and skin.